

SPECIFICATION SHEET

Large Volume Mayer's Hematoxylin

DESCRIPTION:

Mayer's Hematoxylin is a stain with histological applications. This formula is based on the original formulation developed by Mayer and requires slightly longer application periods.

AVAILABILITY:

<u>CATALOG #</u>	<u>VOLUME</u>
TA-060-MH	60ml
TA-125-MH	125ml

PROTOCOL:

1. Rehydrate tissue if needed.
2. Stain nuclei 1-2 minutes with hematoxylin.
3. Place in Bluing Reagent until stain is blue.
4. Rinse in DI water.
5. Coverslip

RESULTS:

Nuclei are blue.

STORAGE:


Store at room temperature.

PRECAUTIONS:

47790 Westinghouse Dr. Fremont CA 94539 USA
Tel: (800) 828-1628 Fax: (510) 991-2826
Email: labvision@labvision.com Website: www.labvision.com



TA-125-MH
Rev 06049
Page 1 Of 1

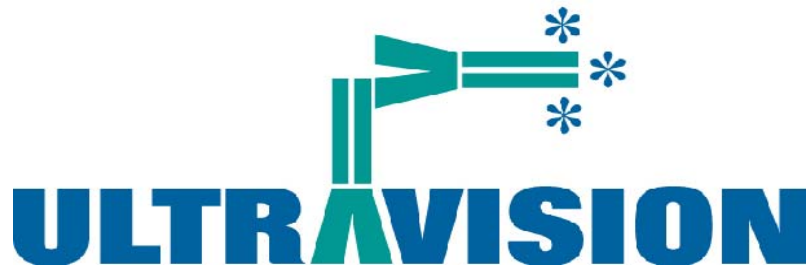

Do not pipet by mouth. Avoid contact with skin and eyes. Wash after use. Observe all federal, state and local environmental regulations regarding disposal.

LIMITATIONS:

This product is available for research use only; not for diagnostic or therapeutic work. Lab Vision Corporation will not be held responsible for patent infringement or other violation that may occur with the use of this product.